


COMPRESSOGRIP PROSTHETIC SHRINKER


Now available with facility logo! See pg. 55!



Fitting Instructions:

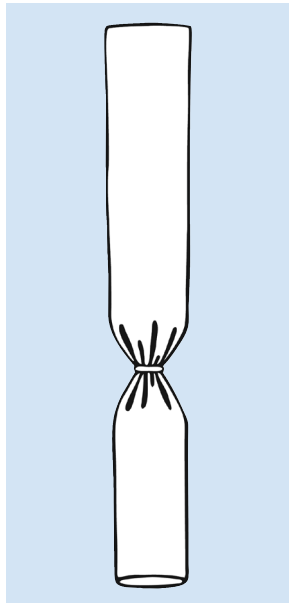

1. Apply the first layer.

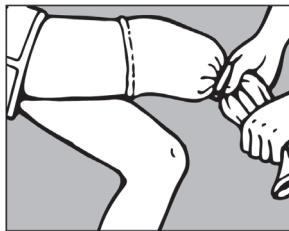

2. Position the ring snugly.



3. Apply the second layer to a point of desired maximum compression.

COMPRESSOGRIP® B/K ELASTIC SHRINKER

The original tubular compression prosthetic shrinker used for control or prevention of edema following amputation surgery, for limb volumetric maintenance, to minimize hypertrophic scarring and to give compressive support for pain relief. Compressogrip Shrinkers assure controlled compression, ease of application and wearer acceptance. Efficacy clinically tested (reprints available).

Compressogrip B/K Shrinker with X-STATIC® – The Silver Fiber® has all the benefits of Knit-Rite's original Compressogrip Shrinker. In addition, the X-STATIC adds permanent antimicrobial benefits to inhibit odor in the garment and effectively wicks moisture away from the skin. All Compressogrip is made from latex-free materials and washable. When determining size, allow for swelling.

White, Beige — 92% polyester, 8% spandex

Polyester/X-STATIC® Grey — 87.4% polyester, 4.6% X-STATIC, 8% spandex



For custom logo shrinkers, see pg.55, or ask customer service or your distributor for program details.

WHITE	BEIGE	GREY/GREY X-STATIC	SIZE	LENGTH	FITS CIRCUMFERENCE
2H0S-02-18	2H0B-02-18	2H-0X02-18	#2	18" (45cm)	5"-7" (13-18cm)
2H0S-02-24	2H0B-02-24	—	#2	24" (60cm)	5"-7" (13-18cm)
2H0S-05-18	2H0B-05-18	2H-0X05-18	#2.5	18" (45cm)	7"-10" (18-25cm)
2H0S-05-24	2H0B-05-24	—	#2.5	24" (60cm)	7"-10" (18-25cm)
2H0S-03-18	2H0B-03-18	2H-0X03-18	#3	18" (45cm)	9"-12" (23-30cm)
2H0S-03-24	2H0B-03-24	2H-0X03-24	#3	24" (60cm)	9"-12" (23-30cm)
2H0S-03-30	2H0B-03-30	—	#3	30" (75cm)	9"-12" (23-30cm)
2H0S-53-18	2H0B-53-18	—	#3.5	18" (45cm)	10"-13" (25-33cm)
2H0S-53-24	2H0B-53-24	—	#3.5	24" (60cm)	10"-13" (25-33cm)
—	—	2H-0XS3-30	#3.5	30" (75cm)	10"-13" (25-33cm)
2H0S-04-18	2H0B-04-18	2H-0X04-18	#4	18" (45cm)	12"-16" (30-40cm)
2H0S-04-24	2H0B-04-24	2H-0X04-24	#4	24" (60cm)	12"-16" (30-40cm)
2H0S-04-30	2H0B-04-30	2H-0X04-30	#4	30" (75cm)	12"-16" (30-40cm)
2H0S-05-18	2H0B-05-18	2H-0X05-18	#5	18" (45cm)	15"-20" (38-50cm)
2H0S-05-24	2H0B-05-24	2H-0X05-24	#5	24" (60cm)	15"-20" (38-50cm)
2H0S-05-30	2H0B-05-30	2H-0X05-30	#5	30" (75cm)	15"-20" (38-50cm)
2H0S-06-18	2H0B-06-18	—	#6	18" (45cm)	18"-24" (45-60cm)
2H0S-06-24	2H0B-06-24	2H-0X06-24	#6	24" (60cm)	18"-24" (45-60cm)
2H0S-06-30	2H0B-06-30	2H-0X06-30	#6	30" (75cm)	18"-24" (45-60cm)

Fits Circumference = Measurement of B/K Circumference 2" Below Medial Tibia Tubercle
Measurement of A/K Circumference 2" From Distal End

NOTE: Compressogrip shrinker fabric gives up to 20-25 mmHg compression at 50% stretch. At 100% stretch the compression is almost doubled. Double layers will almost double compression at both 50% and 100% stretch. Application of more than two layers may exceed the limits of desired therapeutic pressure.

