

COMPRESSOGRIP PROSTHETIC SHRINKER

Now available with facility logo! See pg. 55!

TAPERED COMPRESSOGRIP® B/K ELASTIC SHRINKER

Tapered Compressogrip B/K Shrinkers provide additional width at the top to accommodate fuller thighs and reduce rolling and bunching behind the knee caused by too much compression above the knee. Tubular compression material is used for control or prevention of edema following amputation surgery, for limb volumetric maintenance, to minimize hypertrophic scarring and to give compressive support for pain relief. Compressogrip Shrinkers assure controlled compression, ease of application and wearer acceptance. Compressogrip is made from latex-free materials and washable. When determining size, allow for swelling.

White, Beige — 92% polyester, 8% spandex

Polyester/X-STATIC® Grey — 87.4% polyester, 4.6% X-STATIC, 8% spandex

Patent Pending.

For custom logo shrinkers, see pg.55, or ask customer service or your distributor for program details.

WHITE	BEIGE	GREY/X-STATIC	SIZE	LENGTH	FITS B/K CIRC.
2H05-T4-24	2H0B-T4-24	2H-0XT4-24	#4	24" (60cm)	12"-16" (30-40cm)
2H05-T4-30	2H0B-T4-30	2H-0XT4-30	#4	30" (75cm)	12"-16" (30-40cm)
2H05-T5-24	2H0B-T5-24	2H-0XT5-24	#5	24" (60cm)	15"-20" (38-50cm)
2H05-T5-30	2H0B-T5-30	2H-0XT5-30	#5	30" (75cm)	15"-20" (38-50cm)
2H05-T6-24	2H0B-T6-24	2H-0XT6-24	#6	24" (60cm)	18"-24" (45-60cm)
2H05-T6-30	2H0B-T6-30	2H-0XT6-30	#6	30" (75cm)	18"-24" (45-60cm)


Flat Width	Top
4" (10.2cm)	5" Flat Width; Tapers Last 4" of Length
5" (12.7cm)	6.50" Flat Width; Tapers Last 5" of Length
6" (15.2cm)	7.50" Flat Width; Tapers Last 6" of Length

Fits Circumference = Measurement of B/K Circumference 2" Below Medial Tibia Tubercle


NOTE: Compressogrip shrinker fabric gives up to 20-25 mmHg compression at 50% stretch. At 100% stretch the compression is almost doubled. Double layers will almost double compression at both 50% and 100% stretch. Application of more than two layers may exceed the limits of desired therapeutic pressure.

BONE RINGS


	SIZE
5B705507	7/8"
5B700510	1"
5B700511	1 1/8"
5B7005015	1.5"
5B905506	3/4"

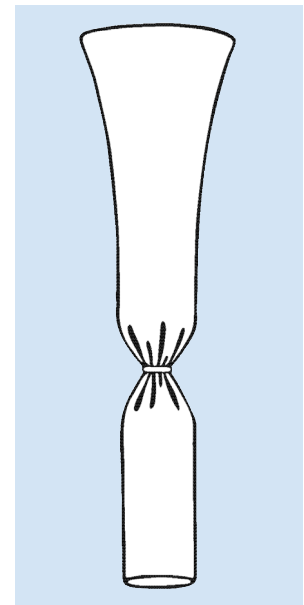

Fitting Instructions:

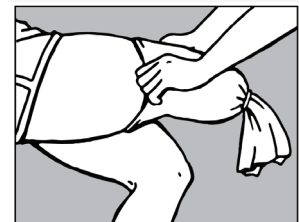

1. Apply the first layer.

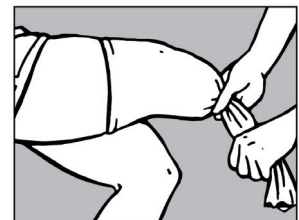

2. Position the ring snugly.



3. Apply the second layer to a point of desired maximum compression.